

# LEPU Non-Contact Forehead Thermometer




## LEPU Non-Contact Forehead Thermometer

The infrared thermometer is intended for the contactless measurement of human body temperature in people of all ages as well as the surface temperature of objects. It may be used by medical professionals or by consumers in a home environment.

### Your benefits at a glance:

- ✓ **SAFE:** Non-contact measurement of body temperature to prevent the spread of bacteria and viruses
- ✓ **FAST:** Temperature is displayed within 1 second
- ✓ **COLOUR INDICATOR:** Coloured backlighting indicates the 3 temperature ranges visually (see pictures on the right)
- ✓ **PRACTICAL:** Ideal for kids and can also be used on unconscious or sleeping people
- ✓ **TWO MODES:** Temperature measurement of human bodies and objects
- ✓ **DATA RECORD:** 99 sets of data can be stored to track changes in body temperature

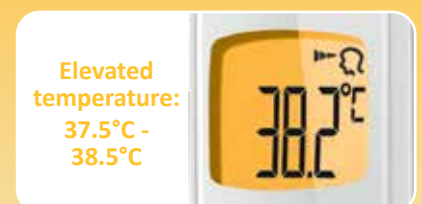

## Product details



Component	Function
1 On/off button	Turns on/off the thermometer. Starts temperature measurement.
2 Memory button	Reviews the stored temperature values.
3 Mode switch button	Switches between body and object mode.
4 Display	Displays measurement value and symbols.
5 Bottom shell	Covers the battery.
6 Sensor	Measures temperature.

## LEPU Non-Contact Forehead Thermometer

By measuring the heat energy emitted from the forehead and calculating body temperature accordingly, accurate readings can be obtained as long as the thermometer is held within 5 cm of the forehead when measuring. Contactless measurement can prevent the spread of bacteria and viruses – safe for children and adults.

Product code	2530104
Name	LEPU Non-Contact Forehead Thermometer
Certification	CE2797
Package contents	Thermometer, instruction manual, warranty card, certificate of conformity, 2 batteries
Application	Contactless temperature measurement within 5 cm of the forehead
Measuring position	Forehead
Measuring distance	0 cm ~ 5 cm
Display	LCD screen
Resolution	0.1°C / 0.1°F
Measuring range	Human body: 32.0°C ~ 43.0°C / object: 0°C ~ 100°C
Accuracy	±0.3°C
Operating temperature	5°C ~ 40°C
Storage temperature	-20°C ~ 55°C
Power supply	2 X 1.5 V AAA batteries
Service life	5 years
Power saving	Auto shutdown when device is idle for 60 seconds

Our **customer service team** is on hand to answer any questions you may have by telephone, via email or in a personal consultation. For details on **your local sales office and regional free phone numbers**, visit [www.nal-vonminden.com](http://www.nal-vonminden.com).

We help you care!

Rapid Tests

Laboratory Diagnostics

Laboratory Service

Consulting & Service

Headquarters

nal von minden GmbH · Carl-Zeiss-Strasse 12 · 47445 Moers · Germany  
 Chief Executive Officers: Sandra von Minden · Roland Meißner · Thomas Zander  
 Commercial reg. Kleve · HRB 5679  
 Tel.: +49 941 290 10-0 · Fax: +49 941 290 10-50

info@nal-vonminden.com  
[www.nal-vonminden.com](http://www.nal-vonminden.com)